

Clinical Changes to RECD Measures: Understanding How the RECD Fits Into the “Big Picture” of Hearing Aid Fitting

Published November 25th, 2015

Danielle Glista, Marianne Hawkins, Sheila Moodie, and Susan Scollie, The National Centre for Audiology, Faculty of Health Sciences, Western University, London, ON

Real ear measurement can allow the clinician to incorporate individual listener’s ear canal acoustics into a hearing aid fitting. Since the introduction of real ear measurement systems, real-ear-to-coupler difference (RECD) measurements became a clinical reality. A hearing aid fitting involves many steps including hearing assessment, hearing aid selection, fitting, verification, and evaluation or outcome measurement. Several recently issued practice standard documents relating to clinical audiology have redefined when and how the RECD should be used to fit hearing aids. Specifically, the American Academy of Audiology (AAA) recommends that when completing pediatric fittings, the RECD be used at the hearing assessment, hearing aid fitting and verification steps¹; recommendation of the RECD measurement as a useful practice in adult hearing aid fitting is made.² RECD measurement is also described as a necessary and routine practice in the prescription of hearing aids for children and adults by the College of Audiologists and Speech-Language Pathologists of Ontario (CASLPO, 2014a, 2014b).^{3,4} Lastly, the new standard on real ear measurement (REM) in North America⁵ (ANSI S3.46 2013) has provided a definition for the RECD measurement. This definition includes the features of the RECD: (1) the coupler type: HA-1; (2) the coupling type: foam tip or earmold; (3) the transducer: high impedance. For the purpose of this article, clinical consideration of coupler and coupling type will be presented in the context of use with a Verifit RE770 RECD transducer. For more detail on issues related to transducer effects and RECD measurement, please refer to Bagatto et al.⁶

The new standardized RECD has recently been incorporated into Verifit hearing aid measurement systems resulting in some small but important changes to how the clinician might set up current hearing aid fittings. Let’s run through some common questions to better understand how the RECD falls into the “big picture” of hearing aid fitting (or on the screen of your hearing aid measurement system). The Audioscan® Verifit system will be used for demonstration; examples will be shown for both the Verifit 1 (VF-1), and the Verifit 2 (VF-2) models.

What is an RECD Used for in Hearing Aid Fitting?

In many fittings, the RECD is used twice before we see the SPLogram and fitting targets on the screen of a measurement system such as the Verifit:

1. To estimate hearing thresholds in sound pressure level (SPL).
HL threshold + RECD + reference equivalent threshold sound pressure level (RETSPL) values = real ear threshold
2. To predict hearing aid output from coupler measures (In Audioscan® systems, the MLE shown

in equation 2 is filtered into the test signal). Coupler SPL or gain + RECD + microphone location effect (MLE) = predicted real ear SPL or gain

Do I Need to Measure an RECD If I Verify Using “Real” Real Ear Measurement?

Although direct verification of the real ear aided response (REAR) does not use the RECD at all, you may still need an RECD to accomplish use #1. The patient’s real ear SPL threshold still gets converted behind the scenes using an RECD in the equation when thresholds are measured with insert phones (the RECD is not applied to thresholds measured with headphones). This may make the displayed thresholds on the SPLogram more accurate for those patients who have unusual ear canal acoustics and result in more accurate target estimation for use in the hearing aid fitting. In comparison, your hearing aid measurement system will need both types of RECDs for coupler based verification to calculate both the real ear SPL threshold and the predicted real ear SPL or gain.

Are RECDs Measured With Foam Tips and Earmolds the Same?

The short answer here is: no. The long answer includes knowing that coupling type (i.e., foam tips or earmolds) does have a direct impact on measured RECD values. This is mainly because the tubing attached to a foam tip is shorter than that attached to most earmolds. RECDs measured with an earmold will often have a high frequency roll-off, when compared to those measured with a foam tip.⁷ For this reason, and because we don’t want to have to complete two RECD measurements, it is best to try and match your coupling selection for the audiogram and RECD. Matching coupling type therefore allows you to complete one RECD measurement that can serve both purposes. Here are three examples of matched protocol for use in hearing aid fitting:

1. Considering babies and young children:

1. Measure hearing thresholds with insert phones coupled to *personal earmold*.
2. Measure RECD once by coupling the RECD transducer to *personal earmolds*.
3. Complete coupler based hearing aid verification.

This first example assumes little or no venting has been included in the fitting. A behind-the-ear (BTE) device would attach directly to an HA-2 coupler (VF-1) or and HA-4 coupler (VF-2) during verification.

2. Considering most adult listeners with a vented/open fitting:

1. Measure hearing thresholds with insert phones coupled to *foam tips*.
2. Measure the RECD once by coupling the RECD transducer to *foam tips*.
3. Complete verification in the real ear.

3. A third example of a matched protocol for use with some adult listeners wearing custom devices (in-the-ear (ITE), in-the-canal (ITC), completely-in-the-canal (CIC) hearing aids) or minimally vented receiver-in-canal (RIC) devices might include:

1. Measure hearing thresholds with insert phones coupled to *foam tips*.
2. Measure the RECD once by coupling the RECD transducer to *foam tips*.
3. Complete coupler based verification (the clinician would need to putty the hearing instrument to the coupler).

What if I Need One RECD Coupling Type, But I Have Only Measured the Other?

For example, what if the audiogram was measured with insert phones coupled to foam tips, the RECD was measured with insert phones coupled to a personal earmold and the device (a BTE, for

example) was verified in the coupler? You might refer to this as a ‘mismatched protocol’, which can result in fitting errors. To reduce such errors, the Verifit applies corrections that happen within the system software. In older versions of Audioscan® Verifit software, the RECD would be used as measured, even if the coupling type was not matched. The current and more accurate approach being used is to automatically generate *foam tip to earmold correction values*.⁷ For example, if the RECD was measured with an earmold coupling, but a foam tip RECD is needed, correction values are applied within the software to generate an individually-predicted RECD for the foam tip coupling. Alternatively, an earmold RECD may be individually-predicted from a measured foam tip RECD by a reversal of this process. To be able to apply such correction values, it is necessary that the clinician define both what type of RECD has been measured, along with some details of the hearing assessment and hearing aid verification setup. If a coupling mismatch is detected in the software, correction values are applied automatically. This brings us to the first important change to be aware of: The VF-1 and VF-2 now require the user to select the RECD coupling methods (e.g., foam tip or earmold) (Figure 1). This is done using the pull-down menu option that best describes the coupling method used. The ‘HL transducer’ menu option is to select the coupling used to measure hearing thresholds, and the ‘RECD coupling’ menu option selects the coupling used during the RECD measurement.

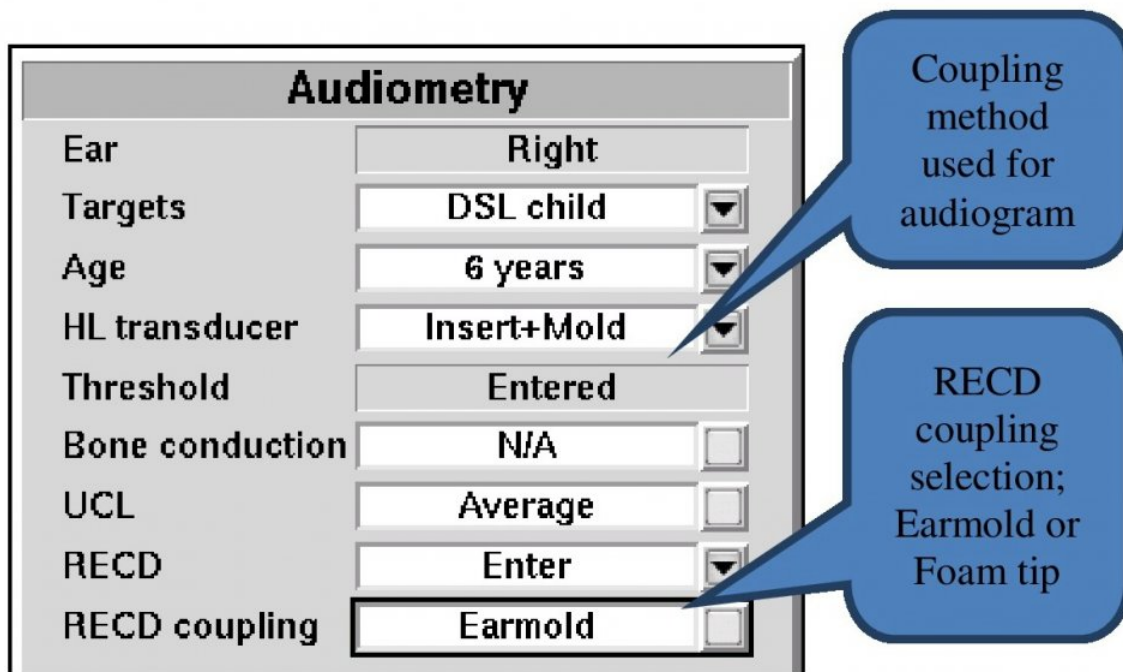


Figure 1.

What Software Versions Include New Earmold To Foam Tip Corrections?

Audioscan® Verifit VF-1 and RM500SL systems (software 3.12 and up).

Audioscan® Axiom (software 1.8 and up).

Audioscan® Verifit VF-2 (software 4.2 and up).

What Coupler Should I Use for Hearing Aid Verification?

Although the new standard includes the HA-1 coupler, coupler type conversions (e.g., between HA-1 and HA-2) are simple, well-understood, and are easily transformed by hearing aid fitting software. The HA-2 is commonly used in verification of BTE devices with the VF-1. This coupler type is no longer used with the VF-2 for hearing aid fitting, except for when completing ANSI testing; a 0.4 cc coupler is used with this system instead since it measures a Wideband RECD

(WRECD) and supports extended-bandwidth verification. Automatic software corrections allow the standard HA-1 RECD to be constructed and reported even in the case when you did not measure your RECD with that coupler type. Therefore, it is possible to complete measurements using different coupler types and have the software reference them back to the “standard” RECD. You will notice that RECDs are labelled as “HA-1” on printed tables of values for this reason.

Application Examples

The following section will go through different verification examples that show possible combinations in coupling and coupler types. These differing combinations can be used when trying to achieve a specific ‘matched’ protocol of RECDs. However, this section also highlights some pros and cons of the different combinations and the need to label each RECD type when setting up a hearing aid fitting. The use of accurate labelling will ensure that the RECD can be corrected and used as needed in the software when there is a mismatch in the protocol. A BTE hearing aid will be used for demonstration purposes. One verification choice is to connect the HA-1 or 0.4 cc coupler to the earmold and BTE device using putty to secure the earmold to the coupler. This procedure requires specific infection control strategies to avoid transferring microbes to the putty, and is not feasible with all types of earmold material. A putty-to-the-coupler example is shown in Figures 2 and 3 for photographs of this, using an HA-1 coupler on the VF-1 (Figure 2) and using 0.4 cc coupler on the VF-2 (Figure 3). The corresponding menu options, illustrating the instrument type and test box coupling method have been included for Figures 2 through 5.

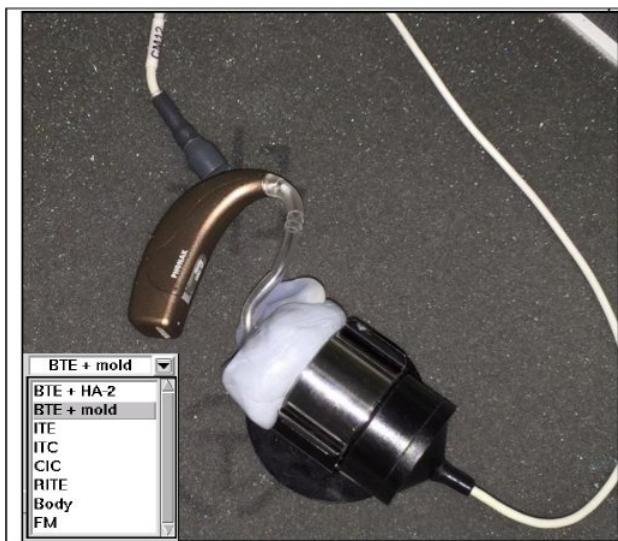


Figure 2.

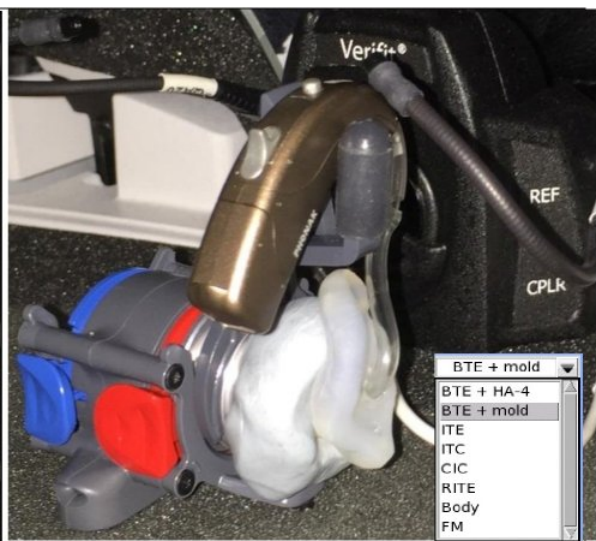


Figure 3.

Alternatively, HA-2 couplers and HA-4 couplers (comprised of the TRIC adaptors and specially designed earmold tubing called the ‘ earmold substitute’) offer the advantage of being easier and cleaner when used with BTE devices because they do not require the use of putty. An example of this method is shown in Figures 4 and 5 for the VF-1 using an HA-2 coupler (Figure 4) and the VF-2 using an HA-4 coupler (Figure 5). If setting up the 0.4cc wideband couplers is new to you, consider watching a tutorial at: <https://www.youtube.com/watch?v=vOAq39Tfb1Y>.



Figure 4.



Figure 5.

What Changes Should I Be Aware of When Measuring or Entering RECDs?

RECD values displayed will depend on the hearing aid system, coupler choice and software version. VF-1 systems require RECD coupler calibration using a 2 cc HA-2 coupler, regardless of software version. The VF-2 uses a 0.4 cc coupler with a TRIC adaptor to connect the RECD transducer. Refer to Figure 6 for an example of the set-up needed to measure the coupler portion of the RECD on the VF-2. Once calibrated, the RECD can be measured on-ear. The VF-1 systems with software prior to version 3.12 will provide HA-2 RECD values. As per the new standard, versions 3.12 and above will convert and display RECDs in HA-1 coupler values when viewing the table of values. VF-2 systems will provide a WRECD which is measured in a 0.4 cc coupler and will convert to HA-1 coupler values when viewing a table of values.

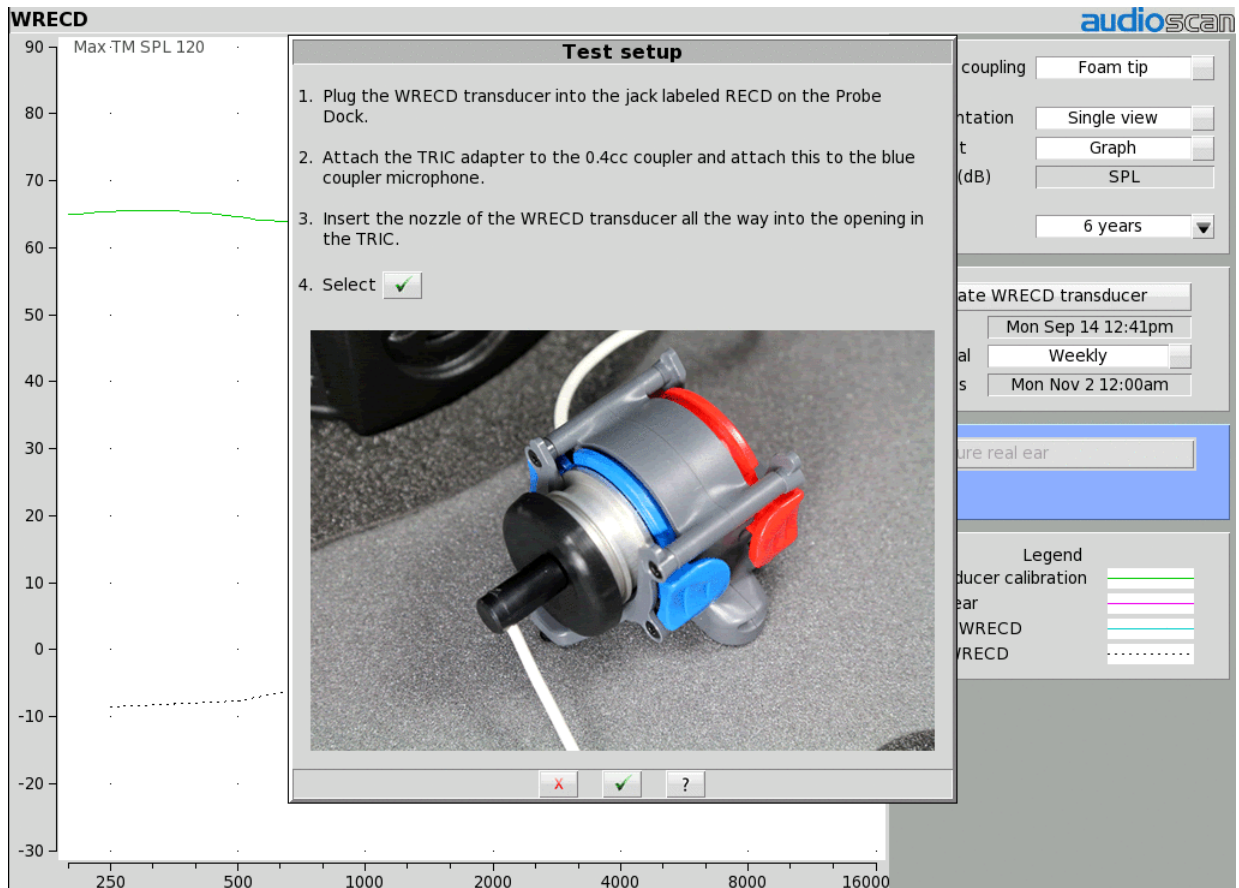


Figure 6.

Previously measured RECDs may be entered into the software, however the system must account for the coupler type used in the original measurement. On-screen prompts appear to assist in the proper selection of a coupler type. Figures 7 and 8 display the new menu options allowing the clinician to enter this information for the VF-1 (Figure 7) and VF-2 (Figure 8). RECD values from older printouts may not be labelled as HA-2 or HA-1. These older RECD were likely measured with an HA-2 coupler and can be entered as such.

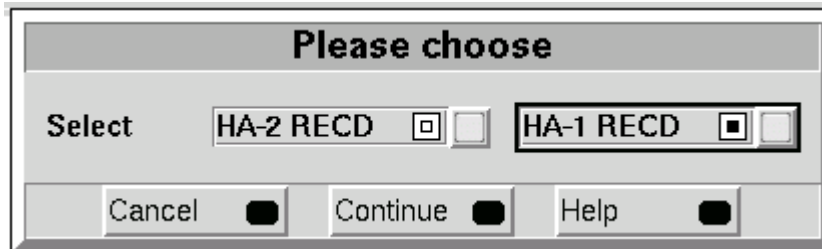


Figure 7.

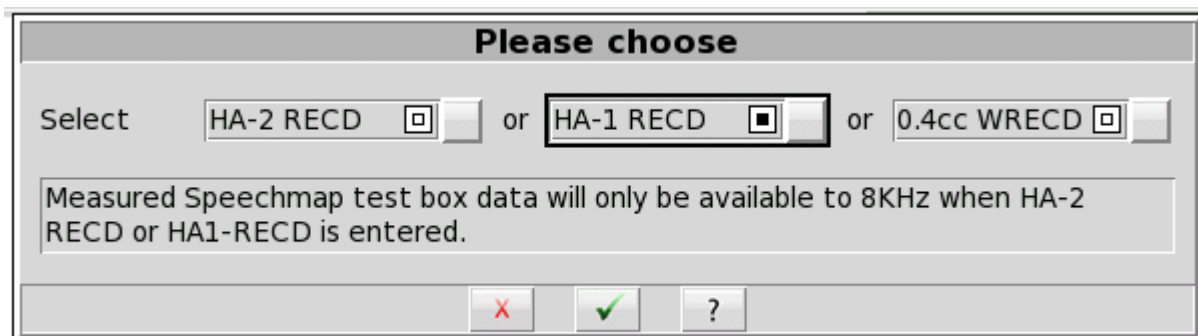


Figure 8.

Where Can I Find A Step-By-Step Guide for Measuring RECDs?

A more detailed description of RECD and WRECD measurements can be seen in this video provided by Audioscan®: <https://youtube/p57vcjTUGYA>. Further details related to specific software updates the VF-2, VF-1 and Axiom can be found at: http://audioscan.com/Verifit2_RECD
http://audioscan.com/Verifit_RECD
http://audioscan.com/Axiom_RECD

Understanding the RECD and how it is used in our hearing aid fitting software allows the clinician to better define hearing aid fitting protocols for individual listeners. Clinicians are faced with many different scenarios that determine how they measure the RECD. For example, children outgrow, lose, and often destroy earmolds, making the measurement of an earmold RECD a challenge. Current hearing aid systems, such as the Audioscan® Verifit, now allow the clinician to label the RECD use and type so a foam tip or earmold RECD can be used to predict the opposite type, and use it as needed.

Acknowledgments

Special thanks to contributions from Chris Stokes-Rees and John Pumford at Audioscan.

References

1. AAA. American Academy of Audiology Clinical Practice Guidelines: Pediatric Amplification. Reston, VA: Author; 2013. Retrieved from http://audiology-web.s3.amazonaws.com/migrated/PediatricAmplificationGuidelines.pdf_539975b3e7e9f1.74471798.pdf.
2. AAA. (2006). Guidelines for the audiologic management of adult hearing impairment. Reston,

- VA Author; 2006. Retrieved from
http://audiology-web.s3.amazonaws.com/migrated/haguidelines.pdf_53994876e92e42.70908344.pdf.
3. CASLPO. Preferred practice guideline for the prescription of hearing aids for children. Toronto: Author; 2014. Retrieved from
http://www.caslpo.com/sites/default/uploads/files/PPG_EN_Prescriptions_Hearing_Aids_Children.pdf
 4. CASLPO. Preferred practice guideline for the prescription of hearing aids to adults. Toronto: Author; 2014. Retrieved from
http://www.caslpo.com/sites/default/uploads/files/PPG_EN_Prescriptions_Hearing_Aids_Adults.pdf.
 5. American National Standards Institute. Methods of measurement of real-ear performance characteristics of hearing aids. Washington; Author; 2013.
 6. Bagatto M, Moodie S, Scollie S, et al. Clinical protocols for hearing instrument fitting in the desired sensation level method. *Trend Amplif* 2005;9(4):199–226.
 7. Moodie S, Pietrobon J, Lindley G, et al. Using the real-ear-to-coupler difference (RECD) within the AAA Pediatric Amplification Guideline: Protocols for applying and predicting earmold RECDs. (in review)